



Product Description


Our high quality Stainless Steel Self Assembly Trolleys are ideal for transporting equipment in cleanrooms, laboratories and other environments that require a focus on hygiene levels

BESPOKE

Do you have a specific furniture design that you cannot find

Utilise our bespoke manufacturing service by contacting us using the details below

FURNITURE

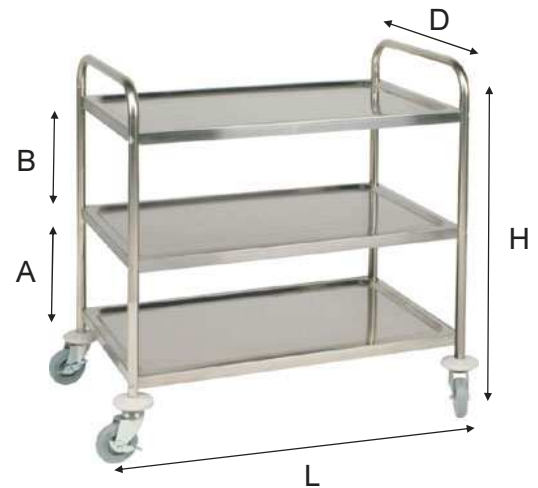

Key Features:

- Manufactured from Grade 201 (Tube) and 430 (Tray) Stainless Steel (See below)
- Carrying capacity of 128kg
- Product comes with a 1 year warranty
- Comes with four quiet swivel castors
- Two lockable castors
- Sound absorbing lining
- Available with 2 or 3 shelves
- Trolley protecting disc bumpers

Delivery Dimensions:

- Product is delivered flat packed
- Dimensions: 710(H) x 145(W) x 430(D) mm
- Weight: 11kg
- Simple self assembly

Size	Shelves	Length (L) (cm)	Height (H) (cm)	Depth (D) (cm)	Distance Between Shelves (cm)	
Small	2	71	81	41	46.5	
Large	2	86	93	54	56.5	
Large	3	86	93	54	A = 29.5	B = 28

Tube: S/S 201 stainless steel:

201 is typical of the 200 series of stainless steels. These steels use manganese to replace nickel in austenitic stainless steels. Low Corrosion - 201 is broadly similar to 304 for corrosion resistance but with some differences in mechanical properties

Tray: S/S 430 0.8mm stainless steel:

Grade 430 is the most common ferritic stainless steel grade in sheet form. It features in many familiar everyday items. Type 430 stainless steel combines good corrosion resistance with good formability and ductility. It is a ferritic, plain Chromium stainless steel with excellent finish quality. Grade 430 also has excellent resistance to nitric attack, which makes it well suited to use in chemical applications.